



Instructions on how to export a CBCT scan from CS 3D Imaging software:



1. Open the scan with the 3D Imaging software.
2. On the side bar under Export, click on “Export the volume”.

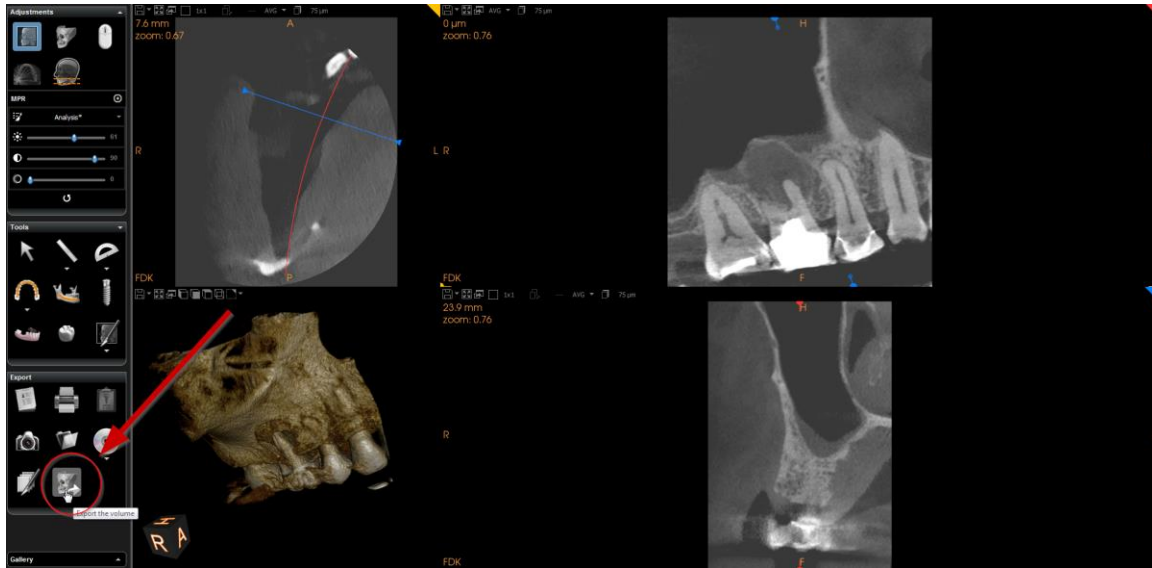

3. Click on “Browse” and create a new folder with the patient name on your desktop: “Patient name”.)

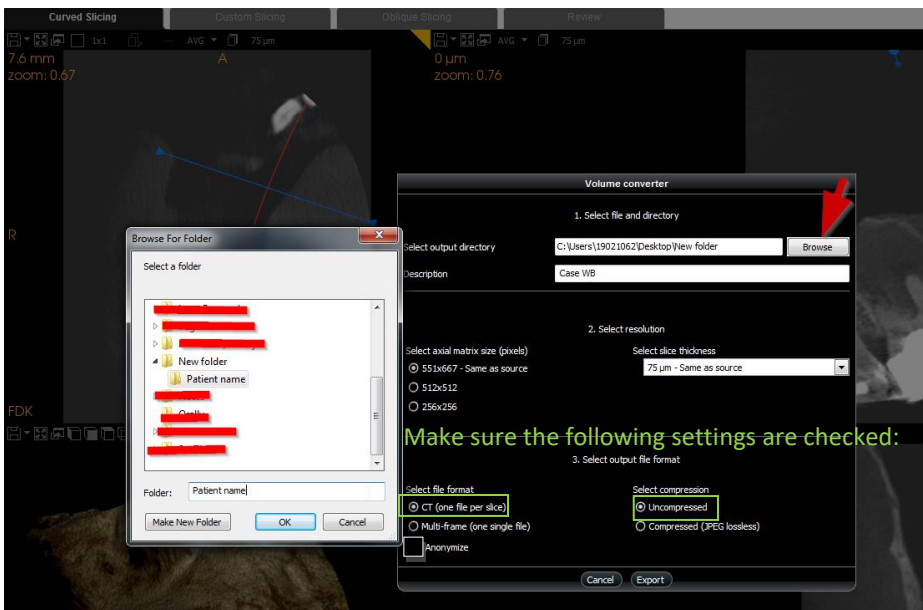

4. Click on Export.
5. Once the export is finished, **right-click** on the folder on your Desktop that has the patient's name on it, go to "**Send to**", and choose "**Compressed (zipped) folder**". This will create a new icon on your Desktop that will have the patient's name on it but the icon will most likely be a folder with a zipper on it. That is the file we need.

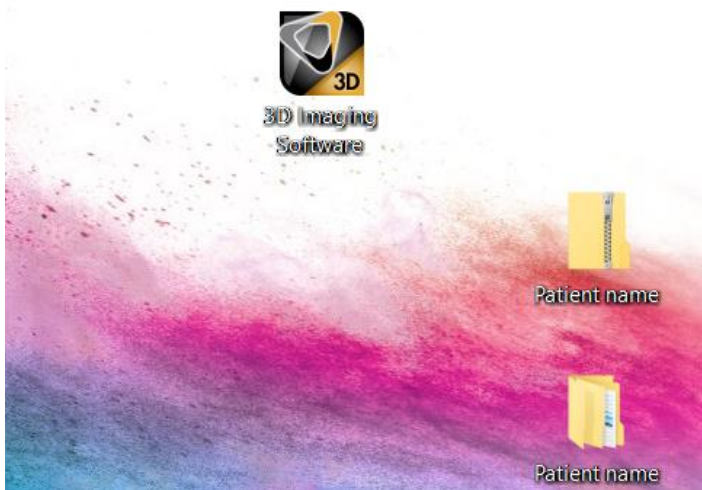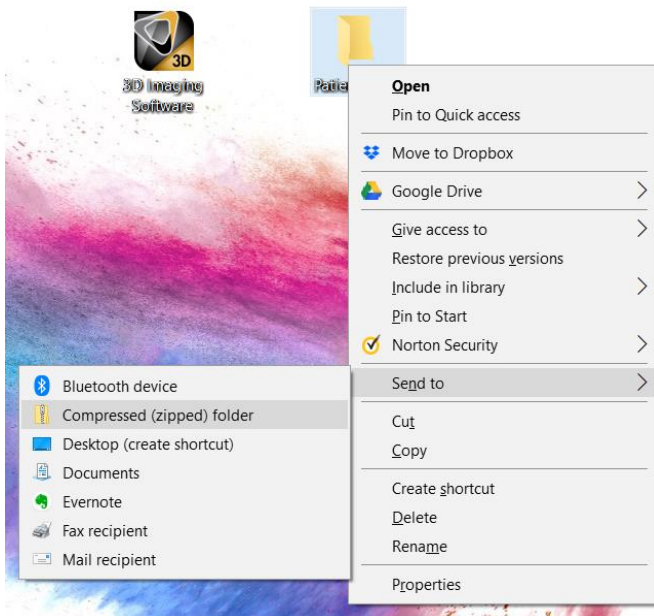

6. Upload the **compressed (zipped) folder** to the HIPAA compliant web portal with the required clinical information.